


AD-I-2071

PIPETTE WASHING STATION

- available for 4x or 8x Autopipettes type AD-K-901
- integration into Positioning Systems possible

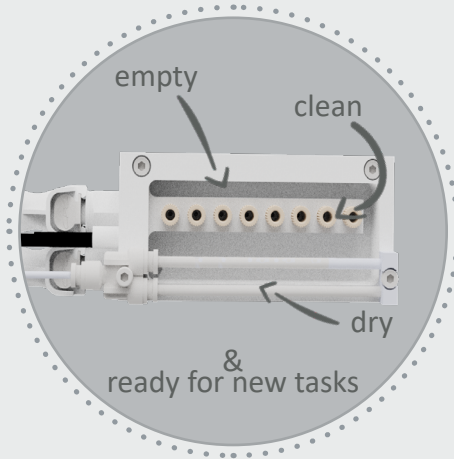

ADVANTAGES


of the built-in Pipette Washing Station:

- Avoidance of Cross-Contamination
- Reduction of Contamination
- Time and Resource Saving
- Extended Lifespan of Pipette Tip
- Improved Result Quality
- Compliance with Quality Standards


See in action:




By utilizing the washing station, you ensure **precision, accuracy**, and an **efficient workflow** in your operations in just three steps:

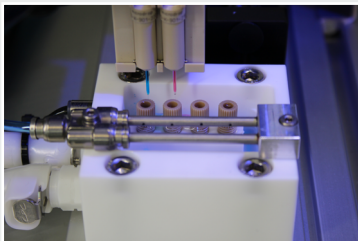

Empty the used liquid(s) after finishing work

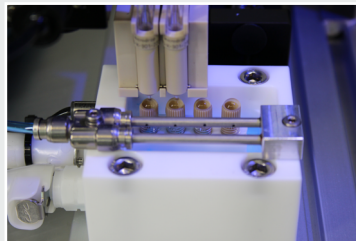

Suck up the cleaning liquid and clean the pipette from inside and outside without cross contamination


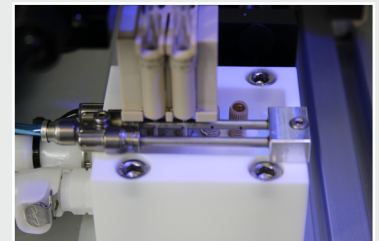


After emptying, the pipette will get blown dry from inside and outside to be ready for new tasks



EASY SETTING

in the included Autodrop Software for one or three different cleaning liquids:

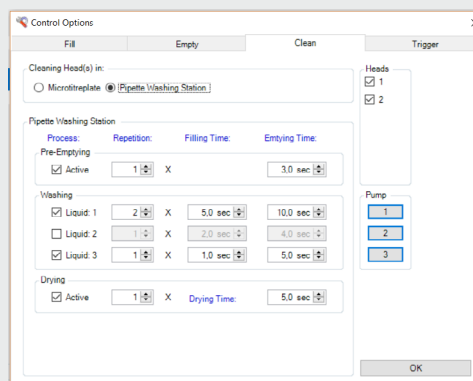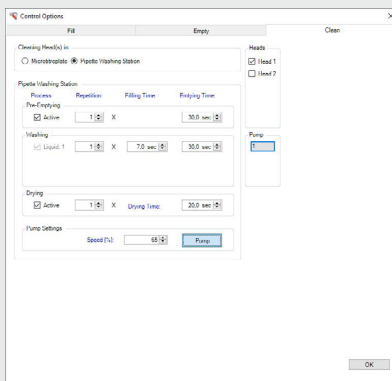

Typical cleaning liquids:

- water
- organic solvent
- IPA
- Ethanol
- Acetone
- Alkaline detergents
- Acid

Intuitive and easy creation of the most suitable cleaning procedure for your application.

We are here for you. Contact us.

microdrop Technologies GmbH
Tycho-Brahe-Kehre 1
22844 Norderstedt
Germany

Tel. +49 40 | 535383-0
Fax +49 40 | 535383-24
info@microdrop.de
www.microdrop.de